


DIA-FLEX FILE



- Made of high-quality Stainless Steel for creating a more fracture-resistant shape and edge
- Enables a complete preparation of even severely curved and narrow canals
- Offers better Cutting Power and Flexibility to shape any root canals easily, safely and very fast
- Available in Sizes 15 through 40 and 21mm, 25mm and 31mm lengths

FLEXIBLE K-FILE

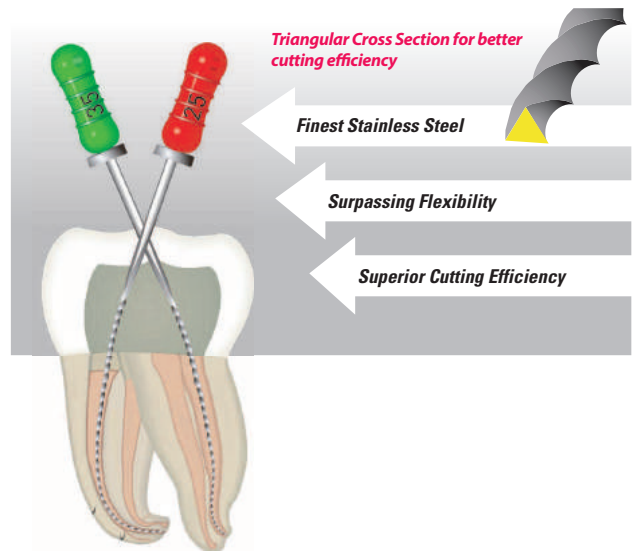

Size	21MM	25MM	31MM
#15	505-103	505-203	505-303
#20	505-104	505-204	505-304
#25	505-105	505-205	505-305
#30	505-106	505-206	505-306
#35	505-107	505-207	505-307
#40	505-108	505-208	505-308
#15/40	505-191	505-291	505-391

(Pkg. of 6)

DIA-SPIRAL FILLER




- Lentulo Paste filler used to evenly distribute root canal sealer and cement in the root canal system
- Conical shape made from grinded wire and made of high grade stainless steel with gradual spirals
- Rather stiff without the safety spiral
- Available in Sizes 25 to 40 and 21mm, 25mm lengths

Size	21MM	25MM
#25	507-105	507-205
#30	507-106	507-206
#35	507-107	507-207
#40	507-108	507-208
#25/40	507-191	507-291

(Pkg. of 4)

DIA-PROBE FILE

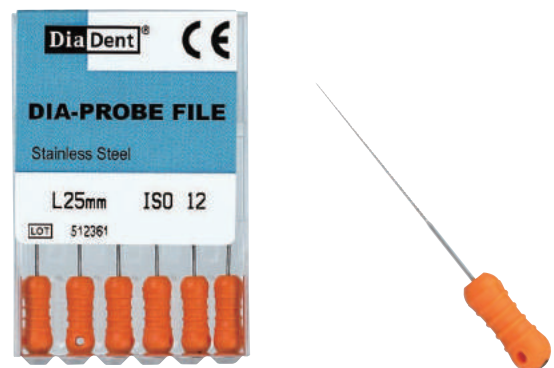

- Canal Probe
- The minimal taper helps to perform intermediate steps in the calcified or curved canals
- Available in Size #12 and 21mm, 25mm lengths

Size	21MM	25MM
#12	506-101	506-201

(Pkg. of 6)

GATE DRILLS

- Used for time-saving preparation of the upper and middle part of non-curved canal or the straight part of curved canal effectively and creating the excellent access for Reamers and Files.
- Enlarge the straight coronal portion of the canal with Gate Drills and then step down with Endodontic Files using the biggest size firstly as the Crown Down Technique.
- Superior Stainless Steel with the improved cutting capacity.
- Circle ring is coded on the shank.
- Recommended Speed: 800~1200r/min

Size	28mm	32mm
#1	531-101	531-201
#2	531-102	531-202
#3	531-103	531-203
#4	531-104	531-204
#5	531-105	531-205
#6	531-106	531-206
#Asst1/6	531-191	531-291

(Pkg. of 6)



PEESO REAMERS

- Used for the preparation of coronal third and picking up the fillings inside the root canal.
- Often used for the post path.
- Excellent stainless steel with a sharper blade to provide the rapid enlargement.
- Circle ring is coded on the shank.
- Recommended speed: 800~1200r/min


Size	32mm
#1	541-201
#2	541-202
#3	541-203
#4	541-204
#5	541-205
#6	541-206
#Asst1/6	541-291

(Pkg. of 6)



ENDODONTIC STAINLESS STEEL HAND FILES

- Root canal cleaning and shaping instruments
- Offers superior cutting efficiency on insertion and withdrawal and procession handling
- Ergonomic color coded handle provides comfort and excellent tactile feedback
- Available in sizes 8 through 80 and 21mm, 25mm & 31mm lengths

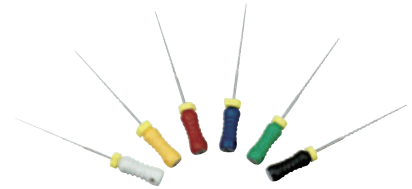

REAMERS



- Made from Stainless Steel wire that is ground to triangular shape and twisted

size	21MM	25MM	31MM
# 8	SS501-101	SS501-201	SS501-301
# 10	SS501-102	SS501-202	SS501-302
# 15	SS501-103	SS501-203	SS501-303
# 20	SS501-104	SS501-204	SS501-304
# 25	SS501-105	SS501-205	SS501-305
# 30	SS501-106	SS501-206	SS501-306
# 35	SS501-107	SS501-207	SS501-307
# 40	SS501-108	SS501-208	SS501-308
# 45	SS501-109	SS501-209	SS501-309
# 50	SS501-110	SS501-210	SS501-310
# 55	SS501-111	SS501-211	SS501-311
# 60	SS501-112	SS501-212	SS501-312
# 70	SS501-113	SS501-213	SS501-313
# 80	SS501-114	SS501-214	SS501-314
# 15/40	SS501-191	SS501-291	SS501-391
# 45/80	SS501-192	SS501-292	SS501-392

(Pkg. of 6)

K-FILES



- Designed to clean and shape root canals
- Produced to give operators a smooth tactile sense inside the canal during instrument

size	21MM	25MM	31MM
# 8	SS502-101	SS502-201	SS502-301
# 10	SS502-102	SS502-202	SS502-302
# 15	SS502-103	SS502-203	SS502-303
# 20	SS502-104	SS502-204	SS502-304
# 25	SS502-105	SS502-205	SS502-305
# 30	SS502-106	SS502-206	SS502-306
# 35	SS502-107	SS502-207	SS502-307
# 40	SS502-108	SS502-208	SS502-308
# 45	SS502-109	SS502-209	SS502-309
# 50	SS502-110	SS502-210	SS502-310
# 55	SS502-111	SS502-211	SS502-311
# 60	SS502-112	SS502-212	SS502-312
# 70	SS502-113	SS502-213	SS502-313
# 80	SS502-114	SS502-214	SS502-314
# 15/40	SS502-191	SS502-291	SS502-391
# 45/80	SS502-192	SS502-292	SS502-392

(Pkg. of 6)

H-FILES



- Ground from Stainless Steel wire
- The flute depth increases with the taper, creating a shaper, more flexible instrument


size	21MM	25MM	31MM
# 8	SS503-101	SS503-201	SS503-301
# 10	SS503-102	SS503-202	SS503-302
# 15	SS503-103	SS503-203	SS503-303
# 20	SS503-104	SS503-204	SS503-304
# 25	SS503-105	SS503-205	SS503-305
# 30	SS503-106	SS503-206	SS503-306
# 35	SS503-107	SS503-207	SS503-307
# 40	SS503-108	SS503-208	SS503-308
# 45	SS503-109	SS503-209	SS503-309
# 50	SS503-110	SS503-210	SS503-310
# 55	SS503-111	SS503-211	SS503-311
# 60	SS503-112	SS503-212	SS503-312
# 70	SS503-113	SS503-213	SS503-313
# 80	SS503-114	SS503-214	SS503-314
# 15/40	SS503-191	SS503-291	SS503-391
# 45/80	SS503-192	SS503-292	SS503-392

(Pkg. of 6)

BARBED BROACHES



- Designed to remove pulp tissue
- Put it straight into the root canal, rotation at 180° and pull out

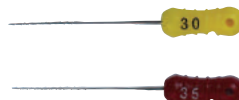

size	21MM	25MM
#20(XXXXF)	504-001	504-201
#25(XXXF)	504-002	504-202
#30(XXF)	504-003	504-203
#35(XF)	504-004	504-204
#40(F)	504-005	504-205
#25/40(XXXF/F)	504-092	504-292

(Pkg. of 6)


SPREADER



- Used for the lateral condensation of root canal fillings.
- Made of Stainless Steel.
- High Resistance to wear and Excellent Cutting Efficiency.
- 6pcs in a blister package

size	21MM	25MM
#15	520-103	520-203
#20	520-104	520-204
#25	520-105	520-205
#30	520-106	520-206
#35	520-107	520-207
#40	520-108	520-208
#15/40	520-191	520-291

PLUGGER

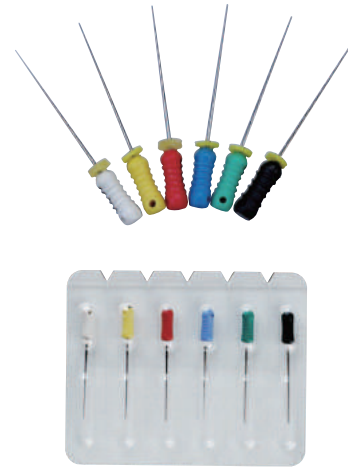

- Used for the vertical condensation of root canal fillings.
- Made of Stainless Steel.
- High Resistance to wear and Excellent Cutting Efficiency.
- 6pcs in a blister package

size	21MM	25MM
#15	521-103	521-203
#20	521-104	521-204
#25	521-105	521-205
#30	521-106	521-206
#35	521-107	521-207
#40	521-108	521-208
#15/40	521-191	521-291

DIA-KONDENSOR

- Having both a nickel-titanium end and a stainless steel end, it can be used in any common endodontic obturation technique including warm and cold condensation.
- Excellent durability.
- Can be supplied 1 piece or a set included each 4 sizes.
- Convenient use with color-coding on the condensers.
- Lightweight and comfort design

size	ITEM No.
#35/70(Green)	551-101
#40/80(Black)	551-102
#50/100(Yellow)	551-103
#60/120(Blue)	551-104
Set (4pcs)	551-191

