

D-LUX

CORDLESS LED CURING LIGHT





D-Lux is a new advanced cordless LED curing light. Made from the highest quality aluminum, D-Lux dissipates heat quickly to prevent overheating problem so that it can produce the maximum light intensity safely and consistently.


Features

- **OLED panel:** unique and advanced one-button control
- A high intensity of **1,600 mW/cm²** offers extremely short curing times of no more than **5 seconds**
- Cordless, highly portable and lightweight (only 155g)
- Unique ergonomic design: provides comfortable grip and easy-to-read OLED panel
- Durable one-piece aluminum design is robust and easy to disinfect
- High-performance lithium-ion battery technology:
 - A new fully charged battery provides 160 x 10-second cures
 - OLED panel conveniently indicates battery life and a countdown timer for the next needed charge
 - Battery can be fully recharged in 90 minutes
- **5 versatile curing modes** for every indications
- Easy to change the battery pack


Product Features

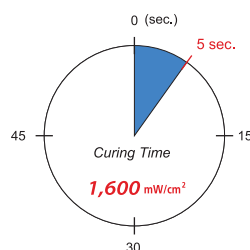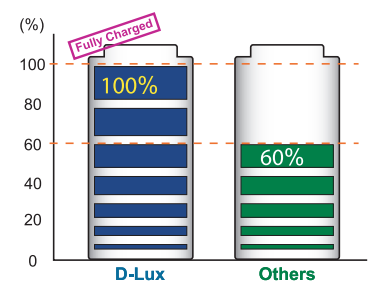Simple & Intuitive
OLED ButtonHigh Intensity
OutputFast Charging
System (90min)


5 Curing Modes

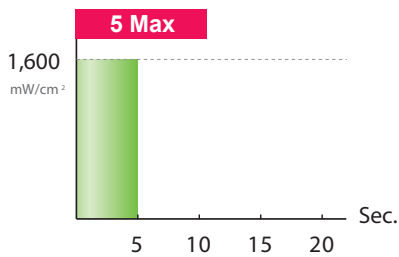


5 Max:
Maximum light intensity output for 5 seconds

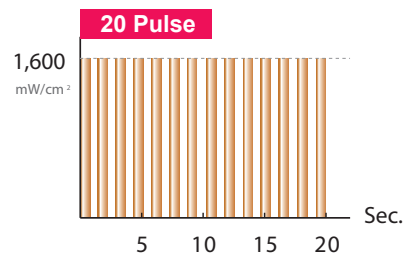


20 Pulse:
Pulses every 0.15 seconds for 20 seconds at maximum light intensity

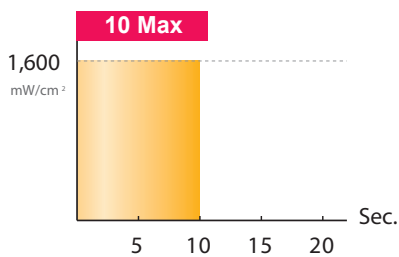

10 Max:
Maximum light intensity output for 10 seconds

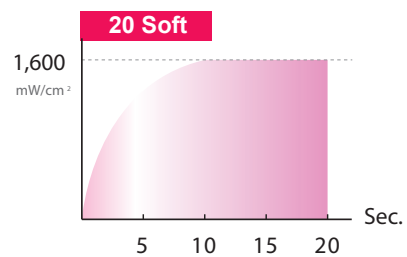

20 Soft:
An initial 10% output of maximum light intensity that gradually increases to full power in 10 seconds, followed by 10 seconds output at full power

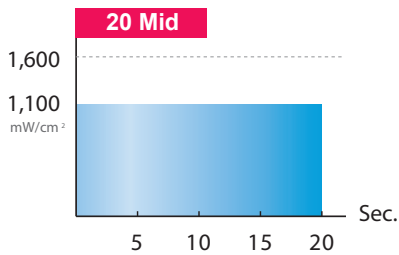

20 Mid:
70% output of maximum light intensity for 20 seconds



Package

- **ITEM #4004-1001 Regular Kit**
 - D-Lux Handpiece & Charger
 - Accessories: Light Guide(11 × 11mm), Plastic Protector, Power Supply and Cord

Device Part Codes

- **ITEM #4004-1110 Handpiece**
- **ITEM #4004-1120 Charger**
- **ITEM #4004-1130 Light Guide(11 × 11mm)**
- **ITEM #4004-1150 Battery Pack**
- **ITEM #4004-1160 Plastic Protector**
- **ITEM #4004-1170 Power Supply (Adapter)**
- **ITEM #4004-1171 Power Cord**